



Theraposture

YOUR INDEPENDENCE IS OUR COMMITMENT

Electrically Adjustable Beds

Freedom Rise

User Manual and Fitting Instructions



Issue 1 - Oct 13



Introduction and General Information.

Thank you for buying the Freedom Rise Bed from Theraposture Ltd. The bed is a high quality product and properly used and looked after, it should give you years of good reliable service. This product is registered as a class 1 medical device with the Medical and Healthcare Regulatory Authority.

In order that you get the best from your bed, please ensure that you and anyone else who is involved with the use of the bed read these instructions carefully.

Please make sure that everyone concerned in using and operating the bed, cleaning it, and cleaning around it, is aware of the section on health and safety. This is particularly important if there are ever children in the area in which the bed is situated.

We would also like to draw your attention to the section on warranty, and particularly the things you can do which may invalidate that warranty.

Please contact us if your requirements change after you have bought the bed.

Please also contact us if you are moving house. In most cases we prefer to move beds that we have supplied. Moving it without contacting us may invalidate your warranty. Enjoy your bed. They make a large contribution to many people's health and well if you are unclear about

anything in this pack, or need some advice, please contact Theraposture on 01373 823030 or email info@theraposture.co.uk. Please always have your sales agreement number (it is in the top right hand corner) and the serial number of the bed (it is on the frame of the bed at the head end) to hand when you call.

Thank you again for buying a Theraposture bed.

If this manual is required in larger print or an alternative format please contact Theraposture 01373 823030 or visit the Theraposture web site at www.theraposture.co.uk

© Theraposture Limited. No part of this document may be transmitted or reproduced, in any form or by any means (system of any nature without written permission of electronic, mechanical, photocopied, recorded or otherwise) or stored in any retrieval system without the prior written permission of Theraposture Limited.

Issue 1 (i) Oct 13 Issue 1 () Oct 13

Contents

Page 1

Intended Use - Limitations of use

Page 2 Product Identification Page

3 Product Description Page 4-5

Health and Safety

Page 6 Product Specification

Page 7 Operating Instructions

Page 8-9 Care and Maintenance

Page 10 Trouble Shooting Page

11 Warranty

Page 12 Disposal Instructions

SECTION 1- FREEDOM RISE

Intended Use

Limitations of Use

The intended use of an electrical bed is to:-

- Assist the user to get into and put of the y

operated

The bed should not be used:-• for which the bed is suitable.To lift or store any object other than a user

- bed.which the bed has been designed.Accommodate only the number of users for •
adult of similar size or a user who may present a risk of entrapmenTo accommodate a
child under 12 or an

- time of use, is within the load capacity of the bed. temperature to suit human habitation. Be
used by a person whose weight at the Be used in a clean, dry environment at a
 - excess of the weight capacity of the bed.has been specificato be able
to understand and operate the hand control. To support a user whose weight is in To support
more than one person unless it By a person who does not have the ability y designed to do so.
-

FREEDOM RISE

- mobility to use the bed effectively and avoid accidents. By a person who does not have sufficient
- damage from the environment Outdoors or where it could be subjected to.

Issue 1

1

Oct 13

Product Identification When making contact with Theraposture about any aspect of the bed,

please have the The identification label can be found at the head end of the bed. Lift the

mattress and it will be located on the metal frame, at fo

FREEDOM RISE

owing information available :-The correct name or model numberThe serial number of the

product If the model number or serial number cannot be found, then the agreement number

will identify the bed. This can be found on the top right hand corner of the sales agreement.

Both can be found on the Identification label, delivery note and warranty certificate. A duplicate identification label to the one fixed to your bed is shown below at Fig 1 below. Make sure the serial number quoted on this page actually relates to the bed in question.

Fig. 1 Identification Label Position



Fix a

Label Here

Serial Number

Issue 1

2

Oct 13

Product Description

Sleep The FreedomRise™ pairs traditional comfort with advanced adjustability for a restful night's sleep. Choose to lay flat similar to the bed you have been sleeping on your entire life or adjust the head and foot sections for a custom sleep position.

Sit The FreedomRise™ can stop in numerous chair positions for the perfect angle to do many of your favorite activities.

FREEDOM RISE

The FreedomRise™ allows you to read, watch TV, play cards, or socialize with friends and family with added comfort.

Rise The FreedomRise™ helps you get up and out of bed to live your life on your own terms. The intuitive two-button operation takes you from a restful night's sleep to a comfortable chair to a standing position to let you go about your daily activities.

Health and Safety

1) The operating instructions for use must

be read before using the bed and must be

revised at a times. 2) The bed must not be used

if faults or defects occur which could result in

risks to people.

3) Before putting the bed into service it, the bed should be checked to ensure that it is functional. These checks should include the following: - y safe and in good working

a) That the mains plug which plugs into the wall is not broken. socket is in good condition

b) The condition of the cable from the mains plug to the control box of the bed is not damaged, broken or trapped in the moving parts of the bed.

c) As smoothly and without noise. movements of the bed operate

d) That there is sufficient clearance around and above the bed to allow the bed to move throughout the full stroke of its movements.
or

4) Electrically powerful motors. The beds are safe if properly used, but care must be taken when using the bed to ensure that damage is not caused when operating these motors. y adjustable beds use extremely

5) in an oxygen-rich environment. (An oxygen mask or nasal bleed may be used) The electric's meet a for this product, BUT the motors are not fully watertight AND are not suitable for use relevant standards

6) Select a suitable location for the bed:

- Free from obstructions.
- With an electrical socket adjacent, so that there is no trip risk from the wire

- If a lifting pole is fitted please ensure that there is sufficient ceiling height available with the bed in the raised position. (Approx 210cm, but it can vary slightly depending on the model of the bed)

7) Whenever the bed is being adjusted,

objects and of the bed's moving parts.

people must be kept clear 8) Ensure that

bedding does not become entrapped within

moving parts of the bed as it is adjusted.

9) The mains cable must be free of the height adjustment mechanism, and be free to move when adjustment is taking place.

10) accessible, but correctly stored off the sleeping platform when adjustments are not being made. Ensure that it is not jammed between any sections of the bed which may cause the buttons to be accidentally operated. Ensure that the hand control is always operated.

11) When the mattress needs replacing, ensure it is replaced with one of exactly the same length, width and thickness.

12) Never drop or insert any object into any opening.

13) Any cords or tubing used on or with this bed MUST be routed and secured properly to ensure that they do NOT become entangled, kinked, squeezed or severed during normal operation of the bed

Issue 1

4

Oct 13

DO NOT use near explosive gases warranty is void. Service is performed on any components by authorized service personnel. If unauthorized persons perform service, the warranty is void. Service is only to be performed by authorized personnel.

14) Cord is plugged into the wall outlet. 15) DO NOT open any actuators, control boxes or pendant cords. 16) DO NOT plug anything into the control box of the bed (i.e. pendant and actuators) while power is on. Possible injury or death may occur due to electrical shock.

Keep a Body weight should be evenly distributed. Possible entrapment hazard may occur if the individual using this bed including mattress and accessories, is over the sleeping surface. Avoid repositioning or transferring out of bed. This includes while assisting the user in you do not use the supplied mattress or if modification are made to the supplied mattress, the bed is being operated, a single segment of the mattress, situations where entire body weight is on to pendant cord being a source for is properly positioned, particularly when disabilities and users with decreased mental acuity should NOT have access to the pendant, entangling user. Children or people with disabilities should not be under the bed at any time. The bed over any power or free of obstructions outlet. **fascinating. However, because of the powerful motors, the bed is not a suitable plaything. Incorrectly used, the bed has entrapment dangers and to the supply wiring.** Children or people with disabilities are **and around the bed, naturally inquisitive, and enjoy playing with Theraposture wishes to draw attention to be in, proximity to the bed. unable to be responsible are in, or likely to be people with disabilities and others** For the above and other reasons: • Specific safety information where children • **the potential dangers of children playing in** such as the one just purchased, moving items. They find beds of this product is 300lbs (136 KG) including accessories being operated. DO NOT permit more than one (1) person side or between parts, especially on the bed at any time. The weight capacity DO NOT let any body parts protrude over the accessories and options. **DO warn visitors of the potential the bed at any time. dangers of the bed. DO ensure that children within the proximity of the bed are well supervised and never left alone. DO NOT allow children to play with the DO warn those with responsibility for the children of the potential dangers of** when the bed

FREEDOM RISE

DO NOT use armrests as handles for moving the bed • DO disconnect the power supply, (when applicable) or, where fitted, disable the hand control using the control box.

Product Specification Actuator Voltage

ITEM Overview

Length (X)

Mattress Width (Z)

Mattress Length (X1)

Mattress Deck Height (Bed) (Y)

Seat Back Angle Range

Seat Bottom Angle Range

Weight of Bed without Headboard or accessories

Maximum Weight Capacity (SWL)

Input Voltage

Features and Accessories

SPECIFICATION

1995 mm (78.5")

915 mm (36") 1905

mm (75")

520mm to 590mm (20.5" to 23.25")

0° to 90°

0° to 19°

FREEDOM RISE

120 Kg (265 lbs)



100-240 VAC, 50/60 Hz

136 Kg (300 lbs)

24 VDC

Standard Features , Foam mattress with Visco topper, High power actuator

Single function handset, Welded frame, Accessories and Options, Fixed Armrests, Prime Mat

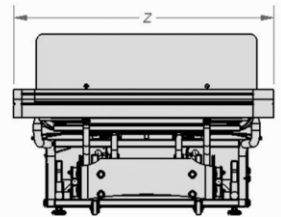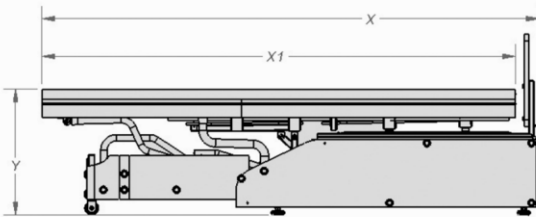

Issue 1

6

Oct 13

Single Function Hand Pendant


••Press top button to convert to chair configuration.Press bottom button to convert to bed configuration.Release button to stop function at any time.

FREEDOM RISE



may be damaged. Secure the power cord to the head end when moving it. the water bed is not in use. Make sure the power cord is not hanging off the product where it **NOTE:** outlet and store the power cord when Disconnect the main power cord from

LEFT/RIGHT REFERENCE GRAPHIC

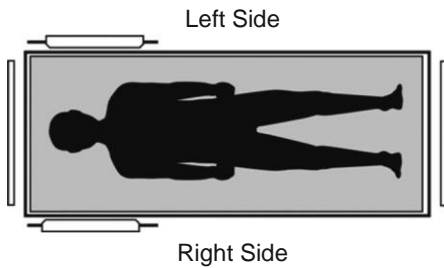

FREEDOM RISE

Cleaning Instructions

- Prior to cleaning unplug power supply cords.
- Make sure all control boxes and pendant are not broken and unplugged. Ensure that NO liquids enter electrical components. housing components are electrical parts (motors,
- DO NOT submerge the bed frame or electrical components.
- Clean the equipment with soap and warm water.
- Solvents, alcohol or petroleum should not be used on the bed.
- Make sure all storing parts are dry before using or

INITIAL INSPECTION

- Inspection of assembled bed components
 - Receipt of
- Check components for obvious damage
- Inspect power supply cords for cuts and/or damage
- Check that actuator cords are connected properly to the controller
- Verify that all are retained by bowtie clips clevis pins are in place and
- Check mattress cover for rips/tears
- Verify proper function

RECOMMENDED INSPECTION CRITERIA

Monthly Inspection

- Inspect all welds for cracks
- Inspect mattress for rips/tears
- Make sure all brackets are tight to the control box/ attachment hardware and

FREEDOM RISE

COMPATIBLE WITH WASHDOWN
TUNNELS **NOTE:** THIS PRODUCT IS NOT

- Check the pendant cord for chafing, cuts or wear and buttons for proper function.

CAUTION

Equipment, property **DAMAGE** or user **INJURY**
can occur during maintenance.

- Check to make sure actuator does not bind at any point throughout its full motion range

INSPECTION

- Inspect wear. Replace or tighten as necessary fasteners for looseness and

Failure to take care of your product

may decrease the life of your product and increase maintenance. Always service the product at the required intervals.

FREEDOM RISE

Issue 1

Oct 13

Care And Maintenance

FREEDOM RISE

Annual Inspection

- Perform a 3 monthly inspections.
 - Check electrical connections for wear or fractures
 - Verify that atight actuator connections are
 - Check actuator cords for chafing, cuts or wear
 - Ensure proper function of any accessories
 - Ensure no tubes are bent
 - Inspect welds on the frame and base assemblies for cracks
 - Verify a fasteners are tight
 - Inspect fasteners for wear or damage
 - Inspect push tube and end connections of actuator for excess wear or bending
 - Verify that aare retained by bowtie clips clevis pins are in place and
-

Troubleshooting Guide

WARNING

Before doing any repairs or maintenance to the product, read aland warnings. The repairs should be done bya skil ed technician. instructions, cautions,

Cause			
Awcotrkuaintorg not	Wcmooar irnnedy eabmceatilogedoonsse	Vilbaoosuacskenelillyghscts oheis wocrk wfrarkyiiinrengg.c.onnCheecctikonsif handfor set	Ranfin daedr n weplf l
Pwrodhileuocpterastatilsng	Tdohewrmn al shut	Bcwifranu ed, itnterms odfe woow ser.ute r .nce w kill s f heor or ctit aok foho sh prmr o ouo r duc t cbs ph we e tt, r ru i od sctiucgioh o hns ot on bed f tim asr e a the ny n	Wa u ap rec th e ttrapl. e 2fgfc c ag theo DO enh
Effect	Possible	Verification	Corrective Action

Warranty Information

1) The Theraposture bed is a quality product, If the bed was purchased new, then it has which should give years of trouble free use. a warranty period as stated on the certificate when the bed was purchased.

~~Damage caused to leads by running them over or catching them in moving parts of the bed.~~
Damage to the mechanism or electrics under the bed by incorrect use of lifting or hoisting equipment. Damage to the upholstery of the bed, however caused.

- 3) This warranty covers parts and labour for faults in the product during that period, unless otherwise stated on the sales agreement. -- Damage caused to a wafurniture. Damage to the electrical supply plug however caused. , or adjacent

- 4) The warranty covers faults which are as the result of a design fault or part failure, whilst being used in accordance with these instructions. 7) The bed is designed to carry a specific weight. This maximum figure is shown on page 6 of this manual. This loading should not be exceeded.

6)-5) Common incidents, which are not covered by warranty, include:-

Accidental

damage to the hand control, for example crushing under a wheelchair or dropping it in

liquid. Incorrect use of the bed, misuse and accidental damage is not covered under the

warranty and there will be no repair. A minor breach of warranty will not affect the continuation of the

warranty. In be a charge for not It is possible to extend the standard warranty,by

one or two years, but only prior to the original warranty expiring. It is possible to

takeadvantage of this offer by contactingTheraposture on 01373 823030 for details.Have

the sales agreement number (it is in thetop right hand corner of the sales agreement)and

the serial number of the bed (it is on the frame of the bed at the head, or page 2 of

this manual) to hand when you identify the product when the call. The above are examples only

and the list is not intended to be exhaustive. This will help to certain

circumstances, however, actions which do fundamental damage to the bed wiof the

warranty period. lead to the invalidation of the remainder

-Liquid spoiling the mattress cover or interior parts, penetrating and degrading the mechanism, electrics, or frame of the bed.

Disposal Instructions - When the product is no longer required, follow these instructions below

- 1) Firstly it should be made safe, so that it is disconnected from the power when it is not in use
- 2) It should also be thoroughly cleaned using disinfectants such as Milton solution to ensure that it is not contaminated.
- 3) If it is to be stored, then please foper the maintenance instructions and then covered in an impermeable dust sheet.ow 1 and 2 above; it should then be greased and oiled as
- 4) If the product has fabric and has been soiled, or if any product has been severely soiled then it should be safely disposed of and not passed on to anyone else.

- 5) If the product is reusable then some thought should be given as to how best to deal with this.
- 6) It is unlikely that Theraposture wisecond hand products so it may be worth contacting the company. buy the bed back, but occasionally have a demand for
- 7) The value of second hand furniture of any sort is very low when compared to the purchase price. Assistive technology furniture is no different and the user should consider any price achieved as a bonus following on from the real reason that it was purchased.
- 8) Theraposture can provide a service inspection + parts if needed.
- 9) Theraposture advise people to advertise locally to see if they have a requirement. Many people recycle such equipment within the family or circle of friends. or on eBay, or perhaps contact their local
- 10) charge. If you want Theraposture to remove the product and dispose of it, there will be a nominal
- 11) involved. Please ring the local authority to discuss this and they may even co_remove the product for a small fee. WEEE regulations but electrical items should be disposed of safely. To dispose of the product, please follow policy as determined by the local authority, particularly where electricals and batteries are charge. Theraposture products are currently excluded from the recycling and environmental ect and

12)If the product is disposed of without them so that the records can be updated to acknowledge this fact.Theraposture involvement, please inform



Theraposture Ltd. Kingdom Avenue, Northacre Industrial Park, Westbury, Wiltshire, BA13 4WE
www.theraposture.co.uk info@theraposture.co.uk Tel: 01373 823030 Fax 01373 826060